



Novel, hand-held, agile surgical operating systems

Nienke Warnaar, Amir Szold

Assia Medical Group, Tel Aviv, Israel

Contributions: (I) Conception and design: All authors; (II) Administrative support: None; (III) Provision of study materials or patients: None; (IV) Collection and assembly of data: All authors; (V) Data analysis and interpretation: None; (VI) Manuscript writing: All authors; (VII) Final approval of manuscript: All authors.

Correspondence to: Amir Szold. Assia Medical Group, Assuta Medical Center, Tel Aviv, Israel. Email: amikisz@gmail.com.

Abstract: Over the past decade, there has been a surge in the development of mechanical and robotized laparoscopic instruments to improve the surgeon's dexterity, precision and ergonomics in a cost-effective manner, as an alternative to the expensive "master-slave" surgical robotic systems currently on the market. Due to the overhead costs, robotic minimally invasive procedures have nil advantage to patients or hospitals when compared to laparoscopic procedures when looking into the cost-effectiveness profile. In addition, to date robotic surgery has statistically not been proven to be superior to conventional laparoscopy for a majority of surgeries. Robotic surgery however continues to grow with increased demand from both patients and surgeons, despite high initial and recurrent costs. As competition expands those costs will inevitably decrease but investing in hand-held robotized laparoscopic devices enabling safe performance of advanced minimally invasive procedures at a much more sustainable cost can enable scaling this technology, especially as every laparoscopic surgeon can easily use those with rapid additional training only to achieve a specific surgical approach. Multiple reports including potential instruments providing those solutions required have been written, however, surprisingly not all of those featured instruments are available today. In this article we review the current status of smart tools for laparoscopic procedures, highlighting new market entries with a potential to change the landscape of minimally invasive surgery.

Keywords: Robotics; minimal invasive surgery; enhanced articulation; hand-held

Received: 03 May 2020; Accepted: 30 October 2020; Published: 20 January 2021.

doi: 10.21037/ales-20-85

View this article at: <http://dx.doi.org/10.21037/ales-20-85>

Conventional laparoscopic instruments are not ergonomic and are restricted to 4 degrees of freedom (compared to 36 in open surgery). This results in severe limitations in performing simple, let alone complex tasks in surgery, holding many surgeons back from engaging in a variety of minimally invasive manoeuvres and procedures. Until recently there were only two categories of laparoscopic instruments: conventional straight manual instruments and a large, console-based, robotic system. Those surgical robotic systems offer increased dexterity, articulation and 3D vision, but at substantial financial costs and logistic complexity, as has been summarized by Perez and Schwaitzberg last year (1).

Even though in a recent study (2) a robotized device did not show superiority compared to conventional laparoscopic instruments in a non-clinical setting, the development

of non-robotic and therefore low-cost laparoscopic instruments that enable better dexterity has recently taken off. Anderson *et al.* reviewed a multitude of mechanical articulating hand-held laparoscopic devices (3). Current articulated mechanical surgical instruments exhibit a wide range of user interfaces, wrist mechanisms and capacities, however, there currently is no clear consensus on what makes an articulated mechanical instrument easy to use. Some articulated mechanical instruments have reached the commercial market and others are under development. As articulated mechanical surgical instruments mature, they have the potential to impact the minimally invasive surgery market by providing some of the capabilities currently only found in robotic systems at a lower cost.

Outside the scope of that review are dexterous

Table 1 Summary of mechanical and robotized hand-held laparoscopic devices

Device	Type	Instrument	DOF	Market availability
FlexDex [®] (FlexDex Inc., Brighton, MI, USA)	Mechanical	Needleholder	6	Worldwide
SILS [®] Hand (Medtronic, Minneapolis, MN, USA)	Mechanical	Interchangeable	7	Worldwide
r2 CURVE (Tuebingen Scientific Medical GmbH, Tübingen, Germany)	Mechanical	Interchangeable	7	Available, mostly Europe
r2 DRIVE (Tuebingen Scientific Medical GmbH, Tübingen, Germany)	Mechanical	Interchangeable	7	Available, mostly Europe
JAiMY [®] (Endocontrol, Grenoble, France)	Robotized	Multifunctional	6	Available, mostly Europe
DEX [™] Robot (Dexterite Surgical, Annecy, France)	Robotized	Interchangeable	7	Available, mostly Europe
HandX [™] (Human Xtensions, Netanya, Israel)	Robotized	Interchangeable	7	Available, mostly Europe
Autonomy LaparoAngle [®] (Cambridge Endoscopic Devices, Framingham, MA, USA)	Mechanical	Needleholder	7	Not available
Realhand [®] (Novare Surgical System, Cupertino, CA, USA)	Mechanical	Interchangeable	7	Not available
Kymerax [®] (Terumo Europe NV, Leuven, Belgium)	Robotized	Interchangeable	5	Not available

DOF, degrees of freedom.

instruments that are partially motorized. A number of hand-held, partially motorized/robotized devices for laparoscopic surgery have been developed, providing additional flexibility in transmitting movement from the user interface to the instrument wrist. These devices however require motors and software, placing them at some cost disadvantage compared to fully mechanical instruments. In *Table 1*, a summary of all available mechanical and robotized laparoscopic instruments can be found. In this short review we will focus on the robotized instruments, as there are the JAiMY[®], the DEX[™] Robot and the HandX[™], in addition to the most widely available mechanical instrument, the FlexDex[®].

The FlexDex[®] is based on a simple and mechanical design, translating the movements of the forearm, wrist and fingers to the tip of the instrument without electrical components (3,4). It provides articulated control and successfully enables suturing in limited spaces. The tool frame is attached as a forearm brace, thus changing instruments may be challenging and time consuming (5-7). In addition, it is compatible with an 8-mm trocar only. Although mechanical, since it is a single use device the cost profile is not low, and currently only a needle holder is available.

The JAiMY[®] has a one-finger control of motorized movements and an ergonomic handle for improved surgical posture, evaluated and proven to be significantly better by

Bensignor *et al.* (8). It is a 5-mm fully articulated reusable instrument, with an unlimited jaw rotation and precise speed control, providing an advantage under ergonomically difficult conditions (9). The device is designed for the purpose of grasping, retracting, mobilizing, dissecting, and suturing of tissues and vessels and is connected by a cable to a control box.

The DEX[™] Robot is a robotized instrument with a grip-type handle, working independently from the shaft, has a complete range of reusable tools and is compatible with any electrosurgical unit. An 8-mm trocar is required. In a small study its use revealed similar results when compared to a conventional instrument for the surgical performance and outcome of an urethrovesical anastomosis. It was therefore concluded that a surgeon's autonomy, dexterity in driving the needle and workload could be improved but with a comprehensive training with the new device. Surgeons did acquire a better body posture using the novel robotic needle holder (10) as were better ergonomics for the surgeon's hand posture observed (11).

In *New Horizons in Laparoscopic Surgery* (2018), a chapter was dedicated to handheld devices (12). Of all the devices that were described to be still in the prototype phase, only the HandX[™] by Human Xtensions from Israel, has been officially launched since. The device received FDA clearance and CE mark. The smart, robotized surgical system integrates all the components required for a modular

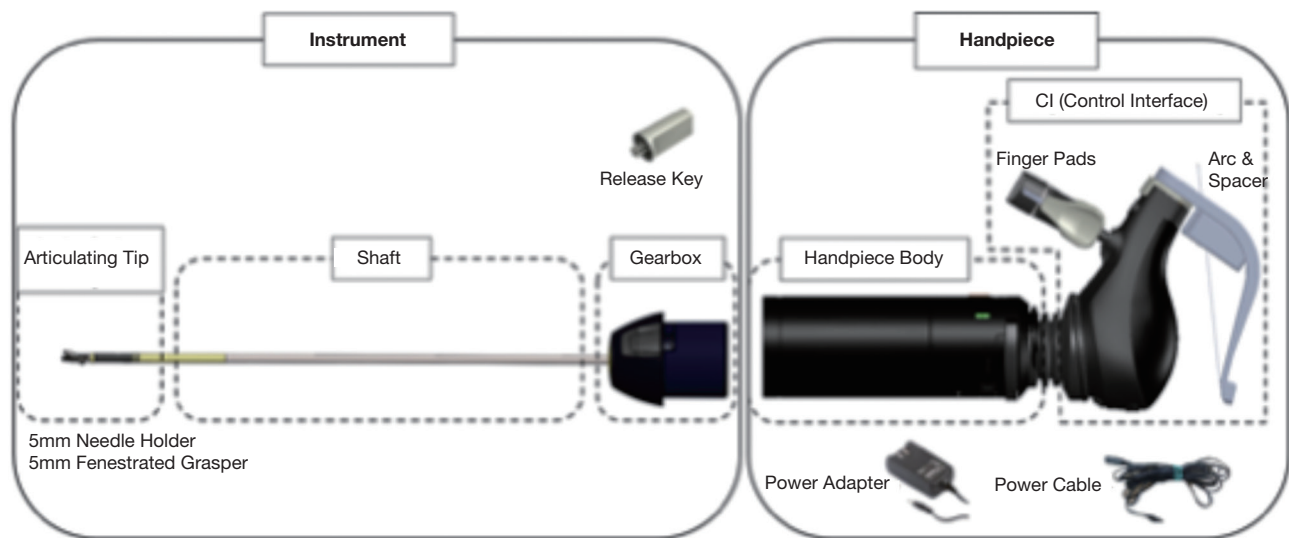


Figure 1 HandX™ device and its component.

platform, of which HandX is the first launched and FDA cleared. It is currently distributed in Europe by Aesculap AG, a subsidiary of B. Braun.

HandX™ (Figure 1) is designed as a light-weight, hand-held device that translates natural unrestricted hand motions into complex movements inside the patient during laparoscopy. The instrument is composed of a sophisticated user interface and a novel, motor driven articulating tool that is controlled by the interface. These components are reusable. The shaft and instrument at the tip are single-patient use and currently a needle holder and grasper are available, although most laparoscopic instruments will be launched soon. The system doesn't require any set up time, and can be easily moved between commercially available laparoscopic 5-mm trocars. Since it was clinically launched over 200 procedures have been performed with the use of this novel device, in multiple countries in Europe as well as in the United States and Israel. The operations included upper gastrointestinal (GI) procedures (sleeve gastrectomies, paraesophageal hernia repairs, gastric bypasses), inguinal and ventral hernia repairs, cholecystectomies, hysterectomies, colectomies, solid organ procedures, thoracic procedures and prostatectomy for benign and malignant disease. The device enabled complex motions and tissue manipulation as well as suturing in difficult angles and in narrow, hard to access spaces. Several clinical trials are now underway to study the use of the device in different settings, and the cost-effectiveness profile of its use.

In conclusion, there are ongoing, interesting and promising developments of smart devices in the area of minimal invasive surgery as an alternative to the currently available robotic systems that are very complicated and costly. In addition, telemanipulated surgical systems lack haptic feedback during the surgical performance, which hand-held devices do provide. Finally, most new motorized instruments are reusable, and the cost-profile of this robot-like dexterity is therefore low. Considering the fact that many of those devices are still in early stages of development, the future for the use of those innovative solutions looks bright.

Acknowledgments

Funding: None.

Footnote

Provenance and Peer Review: This article was commissioned by the Guest Editors (Steven D. Schwartzberg and Rafael Perez) for the series "Advances in Robotic Surgery" published in *Annals of Laparoscopic and Endoscopic Surgery*. The article has undergone external peer review.

Conflicts of Interest: Both authors have completed the ICMJE uniform disclosure form (available at <http://dx.doi.org/10.21037/ales-20-85>). The series "Advances in Robotic Surgery" was commissioned by the editorial office without

any funding or sponsorship. Dr. AS reports equity from Human Xtensions Ltd. during the conduct of the study and is the CMO. NW has no other conflicts of interest to declare.

Ethical Statement: The authors are accountable for all aspects of the work in ensuring that questions related to the accuracy or integrity of any part of the work are appropriately investigated and resolved.

Open Access Statement: This is an Open Access article distributed in accordance with the Creative Commons Attribution-NonCommercial-NoDerivs 4.0 International License (CC BY-NC-ND 4.0), which permits the non-commercial replication and distribution of the article with the strict proviso that no changes or edits are made and the original work is properly cited (including links to both the formal publication through the relevant DOI and the license). See: <https://creativecommons.org/licenses/by-nc-nd/4.0/>.

References

- Perez RE, Schwaitzberg SD. Robotic surgery: finding value in 2019 and beyond. *Ann Laparosc Endosc Surg* 2019;4:51.
- Feng J, Yang K, Zhang Z, et al. Handheld laparoscopic robotized instrument: progress or challenge? *Surg Endosc* 2020;34:719-27.
- Anderson PL, Lathrop RA, Webster RJ III. Robot-like dexterity without computers and motors: a review of hand-held laparoscopic instruments with wrist-like tip articulation. *Expert Rev Med Devices* 2016;13:661-72.
- Awtar S, Trutna TT, Nielsen JM, et al. FlexDex™: A minimally invasive surgical tool with enhanced dexterity and intuitive control. *Journal of Medical Devices* 2010;4:829-39.
- Criss CN, Ralls MW, Johnson KN, et al. A Novel Intuitively Controlled Articulating Instrument for Reoperative Foregut Surgery: A Case Report. *J Laparoendosc Adv Surg Tech A* 2017;27:983-6.
- Criss CN, Jarboe MD, Clafin J, et al. Evaluating a Solely Mechanical Articulating Laparoscopic Device: A Prospective Randomized Crossover Study. *J Laparoendosc Adv Surg Tech A* 2019;29:542-50.
- Gorgen ARH, Araldi M, de Oliveira Paludo A, et al. Laparoscopic pediatric pyeloplasty using the Flexdex® articulating needle driver: step-by-step video. *J Pediatr Urol* 2019;15:421-2.
- Bensignor T, Morel G, Reversat D, et al. Evaluation of the effect of a laparoscopic robotized needle holder on ergonomics and skills. *Surg Endosc* 2016;30:446-54.
- Siri E, Crochet P, Charavil A, et al. Learning Intracorporeal Suture on Pelvitainer Using a Robotized Versus Conventional Needle Holder. *J Surg Res* 2020;251:85-93.
- Sánchez Margallo FM, Sánchez Margallo JA, Skiadopoulou A, et al. Evaluation of the surgical performance, workload and ergonomics during the the urethrovesical anatomises performed with a handheld robotic needle holder. *SAGES*, 2017:abstr 87872.
- Sánchez-Margallo JA, Sánchez-Margallo FM. Initial experience using a robotic-driven laparoscopic needle holder with ergonomic handle: assessment of surgeons' task performance and ergonomics. *Int J Comput Assist Radiol Surg* 2017;12:2069-77.
- Sánchez-Margallo FM, Sánchez-Margallo JA, Szold A. Handheld Devices for Laparoscopic Surgery. In: Ferhatoglu MF. editor. *New Horizons in Laparoscopic Surgery*. IntechOpen, 2018:75-93.

doi: 10.21037/ales-20-85

Cite this article as: Warnaar N, Szold A. Novel, hand-held, agile surgical operating system. *Ann Laparosc Endosc Surg* 2021;6:11.